



Field-proven CAD Technology

Optimal Clinical Flexibility

We understand the constantly changing digital environment and the need for a flexible, compatible CAD solution. R2 ImageChecker delivers proven performance and has been the leading CAD technology for more than a decade.

Using sophisticated software algorithms, ImageChecker searches digital mammograms for potential microcalcifications and masses — characteristics commonly associated with breast cancer. Marks are placed on the regions-of-interest to focus the radiologist's attention. The highly-evolved algorithms, refined to deliver extremely sensitive results without excessive false-positive marks, provide streamlined case review.

Advanced Technology

In addition to detecting regions-of-interest, the ImageChecker algorithms also incorporate anatomic correlation technology. This advanced technology analyzes corresponding findings in the contralateral breast and different views of the same breast.

Flexible Solutions

ImageChecker provides optimal flexibility by allowing you to select sensitivity settings, or operating points, that are most suitable for your needs. ImageChecker ranks findings in order of the prominence of features in a given region. A region will be marked only if the ranking falls above a chosen operating point. Clinical practices can select from three different operating points each for calcification and mass detection, for a total of nine combinations.



ImageChecker supports images from the following digital mammography systems:

Manufacturer	USA	Other
Hologic	Y	Y
GE Medical Systems	Y	Y
Siemens	Y	Y

ImageChecker operates on the Cenova server, providing complete compatibility with DICOM conformant networks.

ImageChecker 9.3 Algorithm Performance¹

The following graphs show plots of algorithm sensitivity vs false-mark rate based on cases with the four standard screening views, with data points for each of the three operating points.

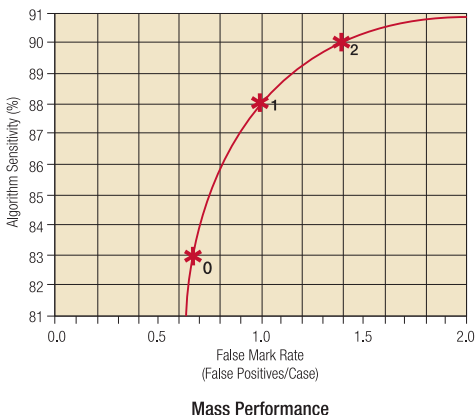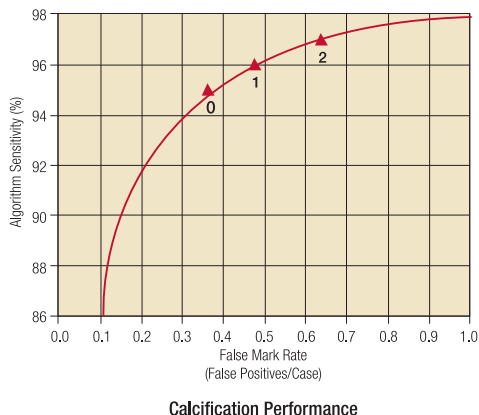

An operating point of "0" corresponds to the lowest sensitivity and the fewest false-marks; an operating point of "2" corresponds to the highest sensitivity and the most false-marks. Practices may select separate operating points for calcifications and masses.

Operating Points

	Calcifications			Masses		
	0	1	2	0	1	2
Sensitivity²	95%	96%	97%	83%	88%	90%
False Marks³	0.36	0.48	0.64	0.68	1.0	1.4
		<i>Default</i>		<i>Default</i>		

Operating Point Combinations

Mass/Calc	0/0	1/1	2/2
Case specificity⁴	48%	35%	24%

Operating Environment

Operating Platform	R2 Cenova
Output Supported	PACS or DICOM workstations

Extended Features

RightOn CAD Marks	Assorted "shaped" markers indicate the types of features that R2 detected
Citra Core ⁵	<ul style="list-style-type: none"> • EmphaSize marks are scaled according to the prominence of features • PeerView Digital provides anatomic outlines of tissue • LesionMetrics provide region specific information, such as lesion <i>Size</i>, <i>Distance to chest wall</i> and <i>Distance to nipple</i>

Performance on Cenova

Maximum Number of Digital Ports	4 (optional)
Case Throughput (four-image)	30-60 cases/hr (nominal) ⁶

¹ Performance derived from analog data. R2 ImageChecker CAD algorithms have been assessed using both digitized film mammograms and digital mammograms. Algorithm performance between the two has been shown to be comparable.

² Data based on 1355 biopsy proven breast cancers (767 mass and 588 calcification cases)

³ Data based on 445 four-view normal cases

⁴ Four film normal cases with no markers

⁵ Available only on advanced SecurView workstations.

⁶ Performance dependent on recommended hardware, network bandwidth and input rate of images

DS-00025 Rev 3 US/International (04/09) © Hologic 2009. All rights reserved. Printed in USA. Specifications are subject to change without prior notice. R2, ImageChecker, RightOn, EmphaSize, PeerView Digital, LesionMetrics, Citra, Selenia and associated logos are trademarks and/or registered trademarks of Hologic, Inc. and/or its subsidiaries in the United States and/or other countries.

HOLOGIC™

The Women's Health Company.

United States / Latin America

35 Crosby Drive,
Bedford, MA 01730-1401 USA
Tel: +1.781.999.7300
Sales: +1.781.999.7453
Fax: +1.781.280.0668
sales@hologic.com
www.hologic.com

Europe

Everest (Cross Point)
Leuvensesteenweg 250A
1800 Vilvoorde, Belgium
Tel: +32.2.711.4680
Fax: +32.2.725.2087

Asia Pacific

Suite 1705, Tins Enterprises Centre
777 Lai Chi Kok Road, Cheung Sha Wan
Kowloon, Hong Kong
Tel: +852.3526.0718
Fax: +1.781.280.0668