SecurMark[®] for eviva[®]

Biopsy Site Marker

Titanium Biopsy Site Marker

Instructions For Use

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Mode d'emploi

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HOLOGIC[®]

SecurMark® Biopsy Site Marker For Eviva® Biopsy Device

Instructions for Use

Please read all information carefully. Failure to properly follow the instructions may lead to unanticipated surgical consequences.

Important: This package insert is designed to provide Instructions for Use for the SecurMark® biopsy site marker for the Eviva® biopsy device. It is not a reference to surgical techniques.

Indications


The SecurMark biopsy site marker is indicated for the permanent radiographic marking of sites in soft tissue.

Contraindications

None known.

Device Description

The SecurMark biopsy site marker is a sterile, single patient use device comprised of a single, biocompatible titanium marker surrounded by a bioabsorbable suture-like material and a deployment device. The deployment device is a hand-held device that delivers the marker from the distal tip. The deployment device consists of a rigid cannula, handle, rigid push rod, plunger, spacer and tip protector. The marker is located at the distal end of the deployment device.

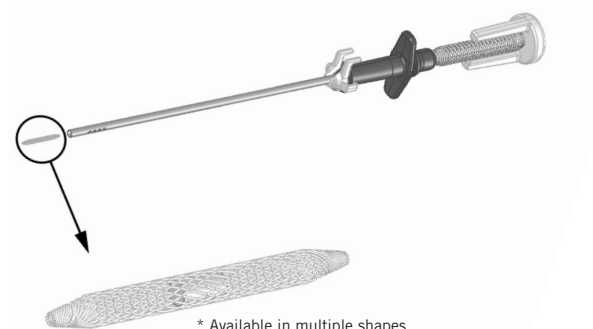Instructions for Use

1. Prior to use of the SecurMark biopsy site marker, inspect the protective packaging and device to verify that neither has been damaged during shipment. If it appears the packaging has been compromised, do not use the device.

2. Remove the SecurMark deployment device from its protective packaging using standard interventional technique.

NOTE: Remove the tip protector prior to use of the device.

NOTE: For use with the Petite aperture biopsy device, remove the spacer from the deployment device by pushing it off with thumb.

3. Remove the biopsy device from the introducer sheath.

4. Place the SecurMark deployment device through the hub of the introducer sheath.

5. Advance the SecurMark deployment device until the handle snaps to the introducer hub.

6. Deploy the marker by advancing the deployment plunger all the way forward until it latches onto the handle.

7. Verify the deployment and proper position of the marker prior to removal of the device with the appropriate imaging modality.

8. Slowly remove the deployment device, or introducer sheath and deployment device, as one unit from the breast and properly dispose.

Warnings and Precautions

- The SecurMark deployment device is not recommended for use within the bore of an MRI magnet.
- The SecurMark biopsy site marker is not recommended for use in patients with breast implants.
- The biopsy site marking procedure should be performed only by persons having adequate training and familiarity with this procedure. Consult medical literature relative to techniques, complications, and hazards prior to performing any minimally invasive procedure.
- The SecurMark biopsy site marker should be used only by physicians trained in open or percutaneous biopsy procedures.
- **Rx ONLY** Caution: Federal (USA) law restricts this device to sale by or on the order of a physician.

- The SecurMark biopsy site marker should be deployed into the cavity created during the biopsy procedure. Deployment into tissue outside of the biopsy cavity is not recommended.
- Marker position relative to established landmarks may change under mammography upon subsequent breast compressions.
- The SecurMark biopsy site marker is not intended to be repositioned or recaptured after deployment.
- Users should take care not to unintentionally deploy the marker.
- Excess hematoma within the biopsy cavity and/or introducer can lead to marker adhesion to the deployment device, increasing the risk of marker drag out.
- Care should be taken to avoid damaging the cannula. Avoid operator or instrument contact with the SecurMark biopsy site marker or the distal end of the deployment device.
- The implanted SecurMark biopsy site marker is magnetic resonance imaging (MRI) conditional. The implanted marker presents no additional risk to the patient or operator from magnetic forces, torque, heating, induced voltages, or movement, but it may affect MRI image quality.
- Minimally invasive instruments and accessories manufactured or distributed by companies not authorized by Hologic, Inc. may not be compatible with the SecurMark biopsy site marker device. Use of such products may lead to unanticipated results and possible injury to the user or patient.
- Instruments or devices which come into contact with bodily fluids may require special disposal handling to prevent biological contamination.
- Dispose of all opened instruments whether used or unused.
- Do not resterilize and/or reuse the SecurMark biopsy site marker device. Resterilization and/or reuse may compromise the integrity of the instrument. This may lead to potential risks of failure of the device to perform as intended and/or cross-contamination associated with using inadequately cleaned and sterilized devices.
- Store the SecurMark biopsy site marker device in a clean and dry area. Avoid storage or handling temperatures above 50°C (122°F).

How Supplied

The SecurMark biopsy site marker device is sterilized by radiation and supplied preloaded for single patient use. Discard into an appropriate container after use.

As Identified on Labels:

QTY

Number of Devices Enclosed.

YYYY-MM-DD

Expiration date is represented by the following:

YYYY represents the year

MM represents the month

DD represents the day

For More Information

For technical support or reorder information in the United States, please contact:

Hologic, Inc.
250 Campus Drive
Marlborough, MA 01752 USA
Phone: 877-371-4372
BreastHealth.Support@hologic.com

International customers, contact your distributor or local Hologic Sales Representative:

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Additional Symbols Used on Labeling

Symbol	Description	Standard
R_x ONLY	Prescription use only	21 CFR 801.109
EC REP	Authorized Representative in the European Community	ISO 15223-1:2016, Reference 5.1.2
CE XXXX	CE Mark with notified body reference number	MDR Regulation (EU) 2017/745
 www.hologic.com/package-inserts	Consult electronic instructions for use	ISO 15223-1:2016, Reference 5.4.3
Translations in Box	Translations in Box	Hologic
	Do not use if package is damaged	ISO 15223-1:2016, Reference 5.2.8
REF	Catalog number	ISO 15223-1:2016, Reference 5.1.6
LOT	Batch code	ISO 15223-1:2016, Reference 5.1.5
QTY	Quantity	Hologic
	Manufacturer	ISO 15223-1:2016, Reference 5.1.1
Made in USA	Made in USA	Hologic
Patents	Patents	Hologic

	Do not re sterilize	ISO 15223-1:2016, Reference 5.2.6
	Do not re-use	ISO 15223-1:2016, Reference 5.4.2
	MR Conditional	ASTM F2503 Reference no. Table 2; 7.4.6.1; Fig 6,7
	Sterilized using irradiation	ISO 15223-1:2016, Reference 5.2.4
	Use-by Date	ISO 15223-1:2016, Reference 5.1.4
	Made in Costa Rica	Hologic
YYYY-MM-DD	Expiration date is represented by the following: YYYY represents the year MM represents the month DD represents the day	Hologic
	Upper limit of temperature	ISO 15223-1:2016, Reference 5.3.6

www.hologic.com/package-inserts

Labeling Translations

SecurMark® for Eviva® Biopsy Site Marker

Product Code: SMARK-EVIVA-10, SMARK-EVIVA-2S-10, SMARK-EVIVA-13, SMARK-EVIVA-2S-13.

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This page is intentionally left blank

This page is intentionally left blank

This page is intentionally left blank

HOLOGIC®

Hologic, Inc.
250 Campus Drive, Marlborough, MA 01752 USA
1-877-371-4372

Hologic BV
Da Vincilaan 5
1930 Zaventem
Belgium
Tel: +32 2 711 46 80

CE
2797