

Glossary of Symbols

Symbol	Standard Reference & Symbol Number	Title of Symbol	Description of Symbol
	EN ISO 15223-1, 5.1.5 ISO 7000, 2492	Batch code	Indicates the manufacturer's batch code so that the batch or lot can be identified.
	EN ISO 15223-1, 5.1.6 ISO 7000, 2493	Catalogue number	Indicates the manufacturer's catalogue number so that the medical device can be identified.
	EN ISO 15223-1, 5.4.4 ISO 7000, 0434A IEC 60601-1, Table D.1, 10	Caution	To indicate that caution is necessary when operating the device or control close to where the symbol is placed, or to indicate that the current situation needs operator awareness or operator action in order to avoid undesirable consequences.
	EN ISO 15223-1, 5.4.3 ISO 7000, 1641 IEC 60601-1, Table D.1, 11	Consult instructions for use	Indicates the need for the user to consult the instructions for use.
	EN ISO 15223-1, 5.2.8 ISO 7000, 2606	Do not use if package is damaged	Indicates a medical device that should not be used if the package has been damaged or opened and that the user should consult the instructions for use for additional information.
	EN ISO 15223-1, 5.3.4 ISO 7000, 0626 ISO 780	Keep dry	Indicates a medical device that needs to be protected from moisture.
	EN ISO 15223-1, 5.1.1 ISO 7000, 3082	Manufacturer	Indicates the medical device manufacturer.
	ISO/DIS 15223-1, 5.7.11 ISO 7000, 6049	Country of manufacture	To identify the country of manufacture of products.
	EN ISO 15223-1, 5.2.7 ISO 7000, 2609	Non-sterile	Indicates a medical device that has not been subjected to a sterilization process.
	ISO 7000, 2794	Packaging unit	To indicate the number of pieces in the package.
	FDA 21 CFR 801.109	Prescription use only	Caution: Federal law restricts this device to sale by or on the order of a physician.
	IEC 60601-1, Table D.1, 20 IEC 60417, 5333	Type BF applied part	To identify a Type BF applied part complying with IEC 60601-1

Warranty

Except as otherwise expressly stated in an agreement between Hologic and its original customer ("Customer"), Hologic equipment ("Equipment") is warranted to the original Customer to perform substantially in accordance with published product specifications for one (1) year starting from the date of shipment ("Warranty Period"). Replacement parts and remanufactured items are warranted for the remainder of the Warranty Period or 90 days from shipment, whichever is longer. Consumable supplies are warranted to conform to published specifications for a period ending on the expiration date shown on their respective packages. Licensed software is warranted to operate in accordance with published specifications. Services are warranted to be supplied in a workman-like manner. Hologic does not warrant that use of Products will be uninterrupted or error-free, or that Products will operate with non-Hologic authorized third-party products.

These warranties do not apply to any item that is: (a) repaired, moved, or altered other than by Hologic authorized service personnel; (b) subjected to physical (including thermal or electrical) abuse, stress, or misuse; (c) stored, maintained, or operated in any manner inconsistent with applicable Hologic specifications or instructions, including Customer's refusal to allow Hologic recommended Software upgrades; or (d) designated as supplied subject to a non-Hologic warranty or on a pre-release or "as-is" basis.

Technical Support and Product Return Information

Contact Hologic Technical Support if the Acessa Handpiece Cable fails to operate as intended. If product is to be returned to Hologic for any reason, Technical Support will issue a Returned Materials Authorization (RMA) number and a biohazard kit if applicable. Return the Acessa Handpiece Cable according to the instructions provided by Technical Support.

If applicable, return used or opened product according to the instructions provided with the Hologic-supplied biohazard kit.

Disposal

Consult your local regulations for disposal and/or electronics recycling as applicable. Do not place into a municipal waste system unless authorized to do so by local authorities.



Rx ONLY

Caution: Federal law (U.S.) restricts this device to sale by or on the order of a physician.

Acessa Handpiece Cable (Model Number 7400) Instructions for Use

Hologic Inc.
250 Campus Drive
Marlborough, MA 01752 USA
www.hologic.com
Telephone : 800-442-9892
<https://www.hologic.com/patent-information>

Hologic, Acessa, Acessa ProVu and associated logos are trademarks and/or registered trademarks of Hologic, Inc. and/or its subsidiaries in the United States and/or other countries.

© 2022-2025 Hologic, Inc.

Product Description:

The Acesa Handpiece Cable provides connectivity of the Acesa Handpiece to the Acesa System Console.

Indications for use:

The Acesa Handpiece Cable is an accessory to the Acesa System Console for use during the Acesa procedure. See the Acesa System User’s Guide (PL-01-0040) for additional information.

Contraindication:

The Acesa Handpiece Cable is contraindicated any time when the use of the Acesa procedure is contraindicated.


Warnings:


- The Acesa Handpiece Cable should be used only by physicians and medical staff who have been trained and have a thorough understanding of the Acesa procedure.
- The Acesa Handpiece Cable is shipped non-sterile

Precautions:

- Do not use if the product is damaged in any way.
- Refer to the Acesa System User’s Guide (PL-01-0040) for additional information on the set-up and use of the Acesa procedure.
- The Acesa Handpiece Cable has not been validated for more than 20 uses.

Preparation and Operation:

1. Remove the Acesa Handpiece Cable from its packaging.
2. Inspect cable for damage. Do not use if it appears to be damaged.
3. CLEAN & STERILIZE THE CABLE PRIOR TO EVERY USE.
4. Ensure that the cable & interconnections are completely dry prior to use.
5. See diagram below: Connect one end to the Acesa Handpiece and the other end to the Acesa Console with this icon: 
6. Ensure that there is a tactile feel when the connector latches or listen to the engaging click of the connector. Do the same for the Acesa Handpiece connector.
7. Follow any other instructions pertinent to the Acesa Handpiece being used.
8. To disconnect the cable, pull straight back on the connector sleeve (see Arrow icon) and then disconnect. DO NOT PULL ON CABLE.

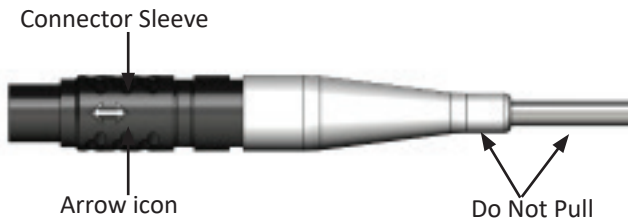

9. Store cable in a dry location.
10. See figure below for connections to the system.



Recommended Cleaning and Sterilization

Hologic, Inc. recommends the following cleaning and sterilization steps as outlined below. Clean and sterilize the Acesa Handpiece Cable before first and every use. Acesa Handpiece Cables must be thoroughly cleaned, rinsed and dried according to the process described in this section before proceeding with sterilization of the device. Thorough cleaning and rinsing are the first and most important steps in the reprocessing of any reusable medical device. Without thorough cleaning and rinsing, it is not possible to achieve effective sterilization of the device.

Cleaning

Rinse contact area of cable thoroughly with running water for two (2) minutes at approximately 43°C (110°F), making sure all foreign matter, blood, mucus, etc. are removed with a soft bristle brush.

Immerse device in Polystica Enzymatic Cleaner, at a concentration of 1/2 ounce cleaner to 1 gallon water, at approximately 43°C (110°F) for three (3) minutes, while wiping and scrubbing cable and connectors with an Endozime Sponge. Follow the detergent manufacturer’s instructions for concentration and other conditions.

Rinse with deionized water for one minute.

Dry with a lint free cloth.

Thoroughly examine all surfaces that have been cleaned and visually inspect the entire device to make sure it is clean.

Table 1. Recommended manual Cleaning Durations

TREATMENT	TIME(MM:SS)	CLEANING SOLUTION
Rinse under running water with a debris removal tool ensuring all soil has been removed	02:00	Water 43°C (110°F)
Immerse device while wiping and scrubbing cable and connector with a debris removal tool.	03:00	Enzymatic Cleaner 43°C (110°F)
Rinse	01:00	Deionized Water

Sterilization

The Acesa Handpiece Cable is provided non-sterile. The cable is sterilized using low-temperature hydrogen peroxide gas plasma technology, selected for its ability to process heat- and moisture-sensitive medical devices quickly, without producing toxic residues or emissions. The cable has been designed to be compatible with ASP STERRAD® 100NX, NX, and 100S systems, as well as STERIS V-PRO® systems.

NOTE: Only one (1) Acesa Handpiece Cable should be run in a STERRAD® NX, 100S, V-PRO® or STERRAD® 100NX load at a time.

The following low-temperature hydrogen peroxide gas plasma sterilization cycles have been validated to result in a SAL of 10⁻⁶ in accordance with applicable standards, including AAMI TIR12:

Table 2. Recommended Sterilization cycles (STERRAD®)

STERILIZATION METHOD	PREPARATION*	STERRAD® STERILIZATION SYSTEM & CYCLE TYPE
Low-temperature Hydrogen Peroxide Gas Plasma Technology	Double wrap utilizing (inner wrap 500 grade, outer wrap 600 grade) FDA cleared Sterilization Wrap	STERRAD® Model 100NX Standard Cycle.
	Place into APTIMAX® Instrument Tray, PC:13837	STERRAD® MODEL NX- Standard Cycle.
	Double wrap with H400 Sterilization Wrap, Halyard Health, PC:68248	STERRAD® Model 100S Short and Long Cycles.

*Device preparation is based on the Acesa validated process. Reference appropriate STERRAD® model instructions for additional information or guidance.

Table 3. Recommended Sterilization Cycles (V-PRO®)

STERILIZATION METHOD	PREPARATION*	STERIS V-PRO® STERILIZER & CYCLE TYPE
Low-temperature Hydrogen Peroxide Gas Plasma Technology	Place into STERIS Sterilization Tray and wrap with H600 OneStep wrap	V-PRO® 1 – Standard Cycle.
		V-PRO® 1 Plus – Lumen and Non-Lumen Cycles
		V-PRO® maX – Lumen, Non-Lumen & Flexible Cycles
		V-PRO® maX 2 – Lumen, Non-Lumen, Flexible and Fast Non-Lumen Cycles
		V-PRO® 60 – Lumen, Non-Lumen and Flexible Cycles
		V-PRO® s2 – Lumen, Non-Lumen, Flexible and Fast Cycles

*Device preparation is based on the Acesa validated process. Reference appropriate V-PRO® model instructions for additional information or guidance.

No toxic residues are left after STERRAD® or V-PRO® sterilization. No aeration is needed for devices sterilized using the STERRAD® or V-PRO® Systems. Specific sterilization questions should be made directly either to Advanced Sterilization Products at: <http://www.asppj.com/>, Or STERIS Corporation at: <http://www.steris.com/>.