

This is a glossary of symbols on Hologic Breast and Skeletal Health products and labels. Use this glossary as a supplement to your product manuals.

Symbol	Description	Standard
	Prescription use only	FDA 21 CFR 801.109
	Not safe for magnetic resonance imaging	ASTM F2503 Reference 6.4.10.1; Fig. 10
	Conditional use for magnetic resonance imaging	ASTM F2503 Reference 6.4.6.1; Fig 6
	Safe for magnetic resonance imaging	ASTM F2503 Reference 7.4.6.1; Fig. 6, 7
	Conditional use for magnetic resonance imaging implant and deployment system	Hologic per ASTM F2503 Reference 7.4.7
	Caution; Radiation	Hologic
	Warning; Radioactive material or ionizing radiation	ISO 7010-W003
	Emergency stop	Hologic
	Provision for the disabled or handicapped persons	ISO 7000-0100
	Footswitch connection	Hologic
	Footswitch cord management plate	Hologic
	Wi-Fi	Hologic

Breast and Skeletal Health Products and Labels

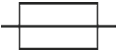
Glossary of Symbols

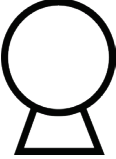
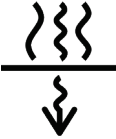
Symbol	Description	Standard
	Danger: Pinch Point	Hologic
	Bluetooth	Hologic
	Battery	Hologic
YYYY-MM-DD	Date format: YYYY represents the year MM represents the month DD represents the day	Hologic
	Circuit breaker	Hologic
Translations in Box	Translations in Box	Hologic
Temp Logger Activated	Temperature Logger activated	Hologic
Temp Logger SN	Temperature Logger Serial Number	Hologic
PN	Part Number	Hologic
REV	Revision	Hologic
Model	Model Number	Hologic
QTY	Quantity	Hologic
SFW SHIPPED	Software version shipped	Hologic
SFW Version	Software version number	Hologic
APPLICATION INSTALLATION	Application Installation	Hologic

Breast and Skeletal Health Products and Labels

Glossary of Symbols

Symbol	Description	Standard
Verified By	Verified By	Hologic
CASE	Case (box or package)	Hologic
VC	Vendor Code	Hologic
Made in USA	Made in United States of America	Hologic
Made in Costa Rica	Made in Costa Rica	Hologic
Assembled in Costa Rica	Assembled in Costa Rica	Hologic
Made in China	Made in China	Hologic
Made in Mexico	Made in Mexico	Hologic
Made in Thailand	Made in Thailand	Hologic
Patents	Patents	Hologic
USA REP	USA Representative	Hologic
LATERAL ARM USE ONLY	Lateral Arm Use Only	Hologic
Q.C. Use Only	Quality Control Use only	Hologic
UNLOAD THIS SIDE	Unload this side	Hologic
RAMP THIS SIDE	Ramp this side	Hologic
MATL	Material	Hologic
	Unlocking	IEC 60417, Reference 5570

Symbol	Description	Standard
	"ON" (power)	IEC 60417, Reference 5007
	"OFF" (power)	IEC 60417, Reference 5008
	Fuse	IEC 60417, Reference 5016
	Protective earth (ground)	IEC 60417, Reference 5019
	Equipotentiality	IEC 60417, Reference 5021
	Alternating Current	IEC 60417, Reference 5032
	Dangerous voltage	IEC 60417, Reference 5036
	Non-ionizing electromagnetic radiation	IEC 60417, Reference 5140
	Type BF applied part	IEC 60417, Reference 5333
	Type B applied part	IEC 60417, Reference 5840
	Computer network	IEC 60417, Reference 5988
	Warning electricity	ISO-7010-W012
	Caution, risk of electric shock	IEC 60417, Reference 6042
	USB	Hologic

Symbol	Description	Standard
	Universal Serial Bus (USB), port/plug	ISO 7000-3650
	Locking, general	IEC 60417, Reference 5569
	Stand-by	IEC 60417-5009
	"ON"/"OFF" (push-push)	IEC 60417-5010
	Earth (ground)	IEC 60417-5017
	Electrostatic sensitive devices	IEC 60417-5134
	"ON" for a part of the equipment	IEC 60417-5264
	"OFF" for a part of the equipment	IEC 60417-5265
	Stand-by or preparatory state for a part of equipment	IEC 60417-5266
	X-ray source assembly	IEC 60417-5338
	X-ray source emitting	IEC 60417-5339
	Radiation filter	IEC 60417-5381

Symbol	Description	Standard
	Do not immerse in any liquid	IEC 60417-5995
IPX8	The equipment of accessory is suitable for continuous immersion in water (up to 1m of submersion for 1 hour)	IEC 60529
 IPX4	Splash-proof	IEC 60529
IPX6	Water projected in powerful jets against the equipment or accessory from any direction shall have no harmful effects	IEC 60529
IP54	Protected against dust limited ingress, no harmful deposits and protected against water splashed from all directions, limited ingress permitted	IEC 60529
IP22	Protected against solid foreign objects of 12,5 mm Ø and greater And Protection against vertically falling water drops when ENCLOSURE tilted up to 15°	IEC 60529
	Medical device	ISO 15223-1:2021, Reference 5.7.7
	Manufacturer	ISO 15223-1:2021, Reference 5.1.1
	Authorized Representative in the European Community	ISO 15223-1:2021, Reference 5.1.2
	Indicates the authorized representative in Switzerland	MedDO; SR 812.213
	Date of manufacture	ISO 15223-1:2021, Reference 5.1.3
	Country of manufacture	ISO 15223-1:2021, Reference 5.1.11

Breast and Skeletal Health Products and Labels

Glossary of Symbols

Symbol	Description	Standard
	Use-by Date	ISO 15223-1:2021, Reference 5.1.4
	Batch code	ISO 15223-1:2021, Reference 5.1.5
	Catalog number	ISO 15223-1:2021, Reference 5.1.6
	Serial number	ISO 15223-1:2021, Reference 5.1.7
	Sterilized using ethylene oxide	ISO 15223-1:2021, Reference 5.2.3
	Sterilized using irradiation	ISO 15223-1:2021, Reference 5.2.4
	Single sterile barrier system	15223- 1:2021 Reference no. 5.2.11 (ISO 7000-3707)
	Double sterile barrier system	ISO 15223-1: 2021 Reference no. 5.2.12 (ISO 7000-3704)
	Single sterile barrier system with protective packaging inside	ISO 15223-1: 2021 Reference no. 5.2.13 (ISO 7000-3708)
	Single sterile barrier system with protective packaging outside	ISO 15223-1: 2021 Reference no. 5.2.14 (ISO 7000-3709)
	Do not re sterilize	ISO 15223-1:2021, Reference 5.2.6
	Non-sterile	ISO 15223-1:2021, Reference 5.2.7
	Do not use if package is damaged	ISO 15223-1:2021, Reference 5.2.8
	Fragile; handle with care	ISO 15223-1:2021, Reference 5.3.1

Breast and Skeletal Health Products and Labels

Glossary of Symbols

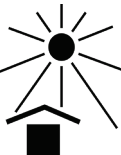
Symbol	Description	Standard
	Keep dry	ISO 15223-1:2021, Reference 5.3.4
	Lower limit of temperature	ISO 15223-1:2021, Reference 5.3.5
	Upper limit of temperature	ISO 15223-1:2021, Reference 5.3.6
	Temperature limit	ISO 15223-1:2021, Reference 5.3.7
	Humidity limitation	ISO 15223-1:2021, Reference 5.3.8
	Atmospheric pressure limitation	ISO 15223-1:2021, Reference 5.3.9
	Do not re-use	ISO 15223-1:2021, Reference 5.4.2
	Consult Instructions for Use	ISO 15223-1:2021, Reference 5.4.3 (ISO 7000 – 1614)
	Non-pyrogenic	ISO 15223-1: 2021, Reference 5.6.3
	Refer to instruction manual/booklet	ISO 7010-M002
	Caution	ISO 15223-1:2021, Reference 5.4.4
	Mass; weight	ISO 7000, Reference 1321A

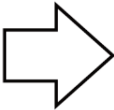
Breast and Skeletal Health Products and Labels

Glossary of Symbols

Symbol	Description	Standard
	This side up	ISO 7000, Reference 0623
	Do not stack	ISO 7000, Reference 2402
	Stacking limit by number	ISO 7000, Reference 2403
	Warning	ISO 7010, Reference W001
	CE Mark European Conformity	MDR Regulation (EU) 2017/745
	CE Mark with Notified Body reference number	MDR Regulation (EU) 2017/745
	Contents	Hologic
	Contains or presence of natural rubber latex	ISO 15223-1:2021, Reference 5.4.5
	Does not contain natural rubber latex	Hologic
	Do not cut	Hologic
	CSA Listed device	Standards Council of Canada (SCC) and the U.S. Occupational Safety and Health Administration (OSHA)

Symbol	Description	Standard
	Medical - General medical equipment as to electrical shock, fire, and mechanical hazards only in accordance with ANSI/AAMI ES 60601-1:2005 (AI:2012), CAN/CSA C22.2 No. 60601-1 (2014), and IEC 60601-1:2012	UL Classification
	ETL Certification	Electrical Testing Laboratories certified
	ETL Certification	Electrical Testing Laboratories certified, Canada and US compliant
	INMETRO & ULBR Mark	UL-BR INMETRO certification (Brazil)
	Recyclable symbol	Unicode, Reference U+267C
	Discard electrical and electronic equipment separately from standard waste. Send decommissioned material to Hologic or contact your service representative.	WEEE Directive 2012/19/EU
	Patient Identification	ISO 15223-1:2021, Reference 5.7.3
	Patient Information Website	ISO 15223-1:2021, Reference 5.7.4
	Health Care Centre or Doctor	ISO 15223-1:2021, Reference 5.7.5

Symbol	Description	Standard
	Date	ISO 15223-1:2021, Reference 5.7.6
	Unique Device Identifier	ISO 15223-1:2021, Reference 5.7.10
	Reject if Dot Black	Hologic
	No pushing	ISO 7010, Reference P017
	No sitting	ISO 7010, Reference P018
	No stepping on surface	ISO 7010, Reference P019
	Warning to use appropriate protection	ISO 7010, Reference W009
	Lifting Hazard	Hologic
	Keep away from sunlight	ISO 7000, Reference 0624
	Contains or presence of phthalates	BS EN 15986, Reference no. A.4

Symbol	Description	Standard
	Contains hazardous substances	ISO15223- 1:2021 Reference no. 5.4.10 (ISO 7000-3723)
	Country code for translation	ISO 3166
	Pinch point	Hologic
	Functional Arrow	ISO 7000, Reference 0251
	Warning: Laser Beam	ISO 7010, Reference W004
	Warning: explosive material	ISO 7010, Reference W002
	Warning: non-ionizing radiation	ISO 7010, Reference W005
	Warning: Floor level obstacle	ISO 7010, Reference W007
	Warning: Toxic material	ISO 7010, Reference W016
	Warning: Crushing	ISO 7010, Reference W019
	Warning: Flammable materials	ISO 7010, Reference W021

Symbol	Description	Standard
	Caution: hot surface	IEC 60417-5041
	Caution Heavy End	Hologic
	Device driver and remote-control connection	Hologic
	Device driver cord management plate	Hologic
	Saline tubing management plate	Hologic
	Option Key	Hologic
	Self-timer	IEC 60417-5879
	Call for maintenance	ISO-7000-0717
	BioVision Camera	Hologic
	Apply specific torque	Hologic
	RCM	AS/NZS 4268, Radio equipment and systems - Short range devices