

Redefining ultrasound access

Introducing the Viera™ portable breast ultrasound





Viera™ portable breast ultrasound **Image accuracy in hand.**

The first wireless, handheld breast ultrasound with
exceptional breast image quality.

Delivers accurate diagnostic images at the point of care.
Ideal for quick diagnostic looks, visual confirmations,
interventional procedures, vascular access and small parts,
**enabling optimization of your clinical workflow and
patient pathway—all at an affordable cost.**



Exceptional imaging

High resolution images transmit wirelessly to your smart device and integrates with PACS.



Streamlined workflow

Handheld, take along anywhere—office, exam room or surgical suite.




Affordable

Add interventional breast ultrasound services for a fraction of the cost of comparable cart systems.

Viera™
Breast ultrasound by **Clarius®**

HOLOGIC®



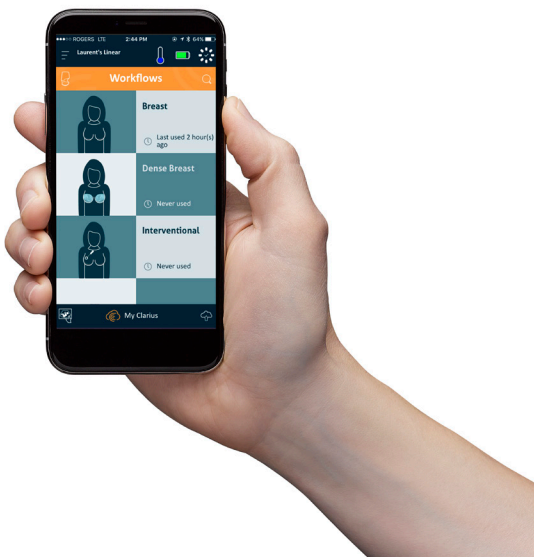

Exceptional Imaging

- High quality images comparable to mid-tier cart systems
- B,M, CFD, PWD and Needle enhancement modes
- 14-4 MHz linear transducer for exceptional resolution down to 7 cm
- 192 elements; 4 software beamformers
- Spatial compounding reduces image noise and speckle



Streamlined Workflow


Breast optimized

- Preset modes for breast, dense breast and interventional procedures simplify workflow
- Customizable annotations package enables efficient documentation
- Breast reporting package for quickly documenting services, ensuring accurate billing and audit



Affordable

- Fraction of the cost of comparable cart systems
- Leverages existing smart device investments
- Reduces scheduling challenges by enabling you to flex any exam space into a procedure room
- Regular and seamless software updates

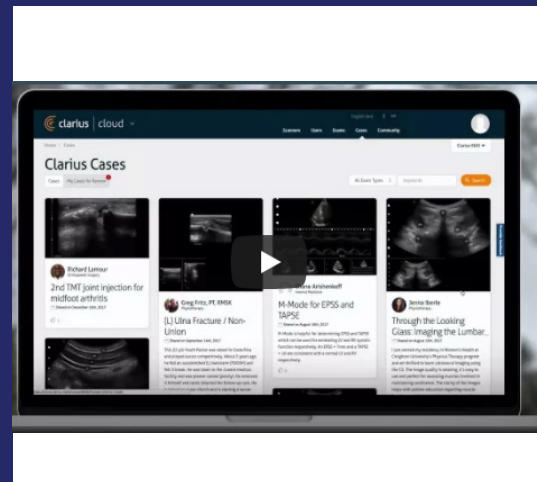

Seamless integration with your existing infrastructure and worry-free maintenance

Viera portable breast ultrasound works with a mobile app that produces high resolution images and transmits them wirelessly to most iOS or Android smart devices.

With many automated imaging controls, including optimal gain and frequency settings, Viera is as easy to use as the camera on your smart phone.



Wireless: Viera uses Wi-Fi Direct to connect to its app on the device (or devices) you choose for imaging. Simply download the free app. Viera does not require Internet access to operate.



Cloud Based: As an authorized user, you can manage your patient's images conveniently and securely with Cloud storage options that integrate with your existing PACS systems or DICOM services of your choice.



Secure: You can rest assured that your patient's data is safe in the cloud, thanks to 128-bit secure SSL encryption and HIPAA-compliant protocols.



Durable: Worry-free maintenance. Built to withstand the demands of clinical settings with a drop safe magnesium shell. No moving parts or cables to break. Sterilizable and submersible up to 3 feet.

“There are no cables, no need to find a plug,
no need for bulky equipment.

I literally walk into a room and start scanning.”

Dr. Kevin Irish
Interventional Radiologist

Streamlined Workflow

The Viera™ ultrasound's exceptional portability and versatility beyond breast applications enable on-demand access to imaging throughout your facility. This will help surgeons, radiologists and breast imagers to efficiently and cost effectively perform a variety of guided interventional procedures such as biopsies, marker placements and wire localizations. In addition, Viera is ideal for quick diagnostic looks, vascular applications, small parts exams and procedures for which you need services at the point of care.



Hologic's interventional products are designed to work together as a complete set of breast biopsy solutions to help you improve efficiency, optimize the patient experience and deliver the best outcomes possible. From image guidance systems to breast biopsy devices, specimen imaging, site markers and accelerated treatment options, Hologic provides the most comprehensive solutions to treat your patients. Be sure to talk with your Hologic representative about which breast biopsy solutions are right for you and your patients.

www.hologic.com | BSHInsideSales@hologic.com | +1.800.442.9892

Hologic Headquarters
United States / Latin America
250 Campus Drive
Marlborough, MA 01752
USA Tel: +1.508.263.2900
Sales: +1.781.999.7453
Fax: +1.781.280.0668
Email: info@hologic.com

Hologic Europe
Everest (Cross Point)
Leuvensesteenweg
250A 1800
Vilvoorde, Belgium
Tel: +32.2.711.4680
Fax: +32.2.725.2087

Hologic Asia Pacific Limited
Unit 01-03A, 13/F,
clar909 Cheung Sha Wan,
No. 909 Cheung Sha Wan Road,
Cheung Sha Wan
Kowloon, Hong Kong
Tel: +852 3748 7700
Fax: +852 3526 0723

Hologic Australia
Suite 402, Level 4
2 Lyon Park Road
Macquarie Park NSW 2113
Australia
Tel: +61.2.9888.8000
Fax: +61.2.9870.7555

PB-00633 Rev.001 (2/18) ©2018 Hologic Inc. All rights reserved. Hologic, Viera, The Science of Sure, and associated logos are trademarks and/or registered trademarks of Hologic, Inc., and/or its subsidiaries in the United States and/or other countries. This information is intended for medical professionals in the U.S. and other markets and is not intended as a product solicitation or promotion where such activities are prohibited. Because Hologic materials are distributed through websites, eBroadcasts and tradeshows, it is not always possible to control where such materials appear. For specific information on what products are available for sale in a particular country, please contact your local Hologic representative.