

The first and only dedicated **prone stereotactic biopsy system with 2D/3D™** imaging.

Introducing the Hologic[®] Affirm[®] Prone Biopsy System



2D/3D™ IMAGING 360° ACCESS



What are you looking for in a prone breast biopsy solution?

The power to conduct both **2D and 3D™ Breast Biopsies** with a single system? **Superior imaging** equal to today's most advanced mammography systems? A **faster, easier way to work? Full 360° access** to the breast? The **power to improve** your patients' biopsy experience?

The Affirm[®] prone breast biopsy system delivers it all.

What's next is here now.

Hologic introduces the world's first and only dedicated prone stereotactic biopsy system capable of both 2D and tomosynthesis-guided procedures. This next-generation solution elevates prone biopsy performance to a level never before possible.



Superior image quality.*

Now you can pinpoint subtle lesions and faint calcifications found during mammography exams. Hologic's Affirm[®] prone system uses the same proven detector technology as our Selenia[®] Dimensions[®] mammography system. It offers a large 14.3 cm x 11.7 cm field of view. And it includes X-ray translucent paddles to deliver high-resolution 2D and 3D™ images.



Streamlined workflow.

Accelerate prone biopsy procedures every step of the way. With automated tube-head motion, programmed biopsy needle parameters, one-click targeting for 3D™ breast biopsies, and more, you can generate high quality images in seconds – for fast procedures and an exceptional patient experience.



Easy 360° breast access.

Quickly access challenging lesions with a fully integrated C-arm that enables 360° breast access to accommodate a wide range of patients. Go from standard to lateral needle approach in seconds to accelerate procedures and ensure you can reach the lesions you're targeting.


* Compared to Hologic MultiCare Platinum system.

Explore the next revolution in prone biopsy.

Identify challenging lesions with a direct-conversion selenium detector and a 14.3 cm x 11.7 cm field of view that is more than 6.5 times larger than existing dedicated prone biopsy systems.

Broaden your view of targeted tissue with X-ray translucent paddles.

Speed up procedures with an automated tube head that sweeps over the breast to acquire all necessary images with a single press of a button.

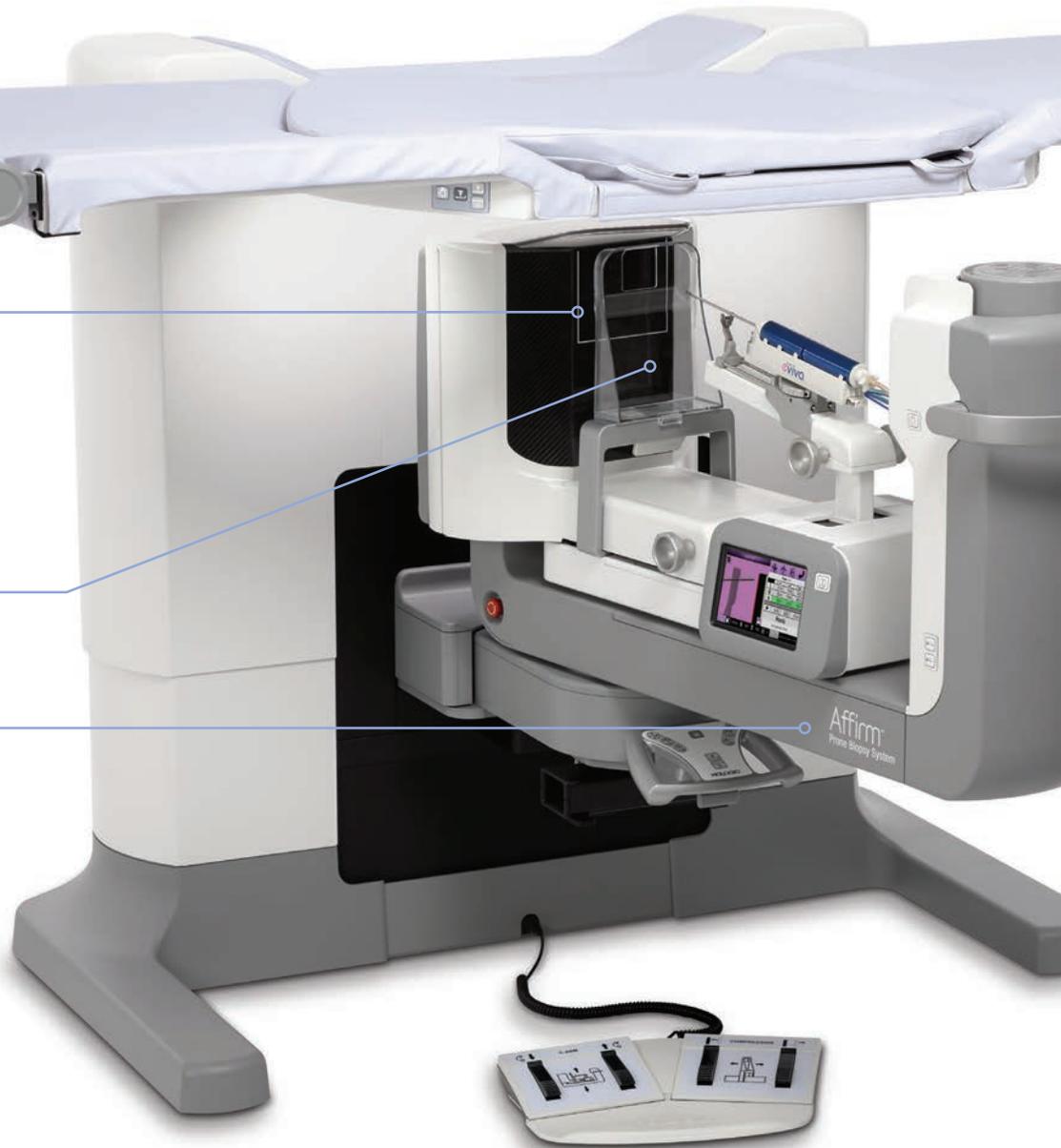

Enhance Patient Comfort.



Maximum comfort package.

Our ergonomically engineered cushions provide custom, targeted pressure-point support for most body types to improve comfort.

Positioning aids.

Arm-through procedure accessories make it fast and easy to comfortably position a patient and access challenging axillary lesions.

Streamlined user experience based on the software found in our Selenia® Dimensions® system.

One-click-targeting when using Hologic 3D™ Breast Biopsy.



Get 360° access to the breast to accommodate most lesion locations – with bidirectional patient positioning and 180° C-arm rotation.

Clearly display the position of the biopsy device to guide the way during procedures via the touchscreen Biopsy Control Module.



Automate calculations – biopsy needle parameters are programmed into the system to display safety margins and relevant distance in real-time.

Innovative lateral needle approach.

Exclusive, fully integrated lateral needle approach that doesn't require additional accessory attachments or taking the patient out of compression.



Hologic's Affirm[®] prone biopsy system delivers an exceptional prone biopsy experience – for everyone.

For radiologists.

- ✔ Use 3D™ imaging to identify and biopsy lesions that only are visible under tomosynthesis imaging.
- ✔ Quickly pinpoint subtle lesions and faint calcifications found during 2D and 3D Mammography™ exams.
- ✔ With a lateral approach and 360° breast access, you can reach challenging lesions and accommodate a wider range of patients; even those with thinly compressed breast.
- ✔ There's no need to remember needle parameters or perform manual calculations, reducing the chance of errors.

For technologists.

- ✔ Quickly identify lesions, making it fast and easy to position patients.
- ✔ Clearly see when a lesion is outside the biopsy window, so you can reposition a patient with confidence.
- ✔ Reduce the number of manual steps needed to complete a biopsy.
- ✔ Perform QC just once a day instead of before every procedure to boost productivity.
- ✔ Programmed needle parameters eliminate manual calculations and reduce the chance of mistakes.



For patients.

- ✓ The faster you can identify, target and access lesions for biopsy, the less time under compression for your patients. The Affirm prone biopsy system helps improve the speed and efficiency of many prone biopsy procedures.
- ✓ According to a recent Hologic study, doctors prefer prone breast biopsy and believe it's more comfortable for a majority of patients.*

Which Affirm system is right for your patients? How about both?

The Hologic Affirm system portfolio now includes both upright and prone breast biopsy guidance systems, so you can provide minimally invasive procedures to the largest group of patients – without compromise.

Our breast biopsy systems share many of the same advanced technologies, software, features, and operations to deliver superior imaging, streamlined workflow, and maximal access. It's the reason so many doctors choose Hologic systems to help improve patient comfort and satisfaction.

Be sure to talk with your Hologic representative about which breast biopsy solutions are right for you and your patients.

* Data on file 2015. Compared to the Hologic MultiCare® Platinum system.

Hologic is a global company using **The Science of Sure** to improve lives by innovating and bringing to market the best diagnostic, detection and surgical technology in the world. We integrate The Science of Sure into everything we do to help healthcare professionals discover, diagnose, and treat their patients with progressive certainty to deliver ever-greater **peace of mind**.

www.hologic.com | info@hologic.com | +1.781.999.7300

**Hologic Headquarters
United States / Latin America**
250 Campus Drive
Marlborough, MA 01752
USA Tel: +1.508.263.2900
Sales: +1.781.999.7453
Fax: +1.781.280.0668
Email: info@hologic.com

Hologic Europe
The Corporate Village
Da Vincilaan 5
Building Caprese,
8th floor
1935 Zaventem
Belgium
Tel: +32.2.711.4680
Fax: +32.2.725.2087

Hologic Asia Pacific
7th Floor, Biotech Centre 2
No. 11 Science Park West Avenue
Hong Kong Science Park
Shatin, New Territories, Hong Kong
Tel: +852.3748.7700
Fax: +852.3526.0723

Hologic Australia
Suite 402, Level 4
2 Lyon Park Road
Macquarie Park NSW 2113
Australia
Tel: +61.2.9888.8000
Fax: +61.2.9870.7555

PB-00440 Rev.004 (8/17) Hologic Inc. ©2017 All rights reserved. Hologic, 3D, 3D Mammography, Affirm, Dimensions, Eviva, Selenia, The Science of Sure, and associated logos are trademarks and/or registered trademarks of Hologic, Inc., and/or its subsidiaries in the United States and/or other countries. This information is intended for medical professionals in the U.S. and other markets and is not intended as a product solicitation or promotion where such activities are prohibited. Because Hologic materials are distributed through websites, eBroadcasts and tradeshow, it is not always possible to control where such materials appear. For specific information on what products are available for sale in a particular country, please contact your local Hologic representative.